


Stress Test ECG




Stress Test ECG "CardioScan Stress Test"

CardioScan Stress Test Version 4.0 = High Resolution ECG+Superior Screen Displays

- Higher Sampling Rate (4096 Hz) and Superior Resolution (12bits) for best ever ECG signal quality
- Custom ECG Sampling Screen Displays for clear and high quality viewing of real time ECG, Superimposed Reference/actual Averaged ECG Complex, Full Disclosure, in 1, 3, 6 or 12-lead mode
- Instant Real Time ECG print outs or ECG Event Marker for later ECG review
- DM Software Holter algorithm for Heart Rate Triggering, ECG Averaging, wandering baseline correction, antialiasing and noise filtering technique
- User settings for ST analysis from J-point or ST point, Absolute or Delta ST Method, ST elevation and depression criteria
- Actual ECG versus Reference ECG or prior stage ECG comparison during test
- Colour coded Arrhythmia analysis from DM Software Holter ECG technology
- Auto, Hold, Manual modes for Treadmill and Ergometer control
- Auto or Manual NIBP measurement input
- Full Disclosure ECG Review, Heart Rate Recovery, individual ST Graph with ECG replay with superimposition, multi trends and graphs, ST reanalysis
- Multiple well known Protocols integrated and unlimited creation of new protocols
- Multi Ergometers, Treadmills, VO2 interfaces, 2 Analog ECG outputs
- Windows 2000, NT, XP, Vista, Tablet PC Compatible
- E-mail, PDF, ASCII, GDT, HL7, Dicom interfaces...

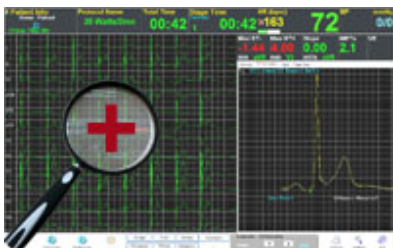

Stress Test Realtime Screen Display

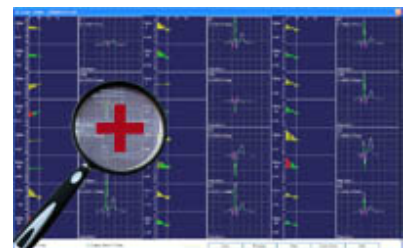

ST Analysis Review Post Test

Hardware



12-Lead Blue-tooth Transmitter: Front & Back



Blue-Tooth Receiver



Patient Cable

Specifications and Features table

Model	12-Lead ECG Blue Tooth Module	12-Lead ECG Blue Tooth Transmitter
Leads	12-Lead, 10 electrodes	12-Lead, 10 electrodes
Sampling Rate	12 Bits, 4096Hz	12 Bits, 4096Hz
Dimension	Trans. 13.5 x 7.5 x 2.4 cm Receiv. 13.9 x 12 x 3 cm	13.5 x 7.5 x 2.4 cm
Weight	Emet. 110 g récept. 240 g	110 g
SAECG * & VCG *	Yes	Yes
FCG *	Yes	Yes
Power	Trans. 4 x AA batt. (6.0V) Receiv. USB Port (5.0V)	4 x AA batteries (6.0V)
ECG Download	Blue Tooth Transmission 2,4 GHz 100 M	Blue Tooth Transmission 2,4 GHz 100 M
Compatibility	CardioScan Resting 12-Lead CardioScan Stress Test	CardioScan Resting 12-Lead CardioScan Stress Test

High Resolution ECG sampling devices for Desktops, Laptops, Tablet PC.

■ New Generation 10m range Blue Tooth Transmitter/Receiver Module for wireless ECG ensuring better ECG quality, Ideal Patient Isolation and Freedom of Movement and Intervention, for use with Desktop, Laptop or Tablet PC!

■ Mini USB module for direct wire link to Desktop, Laptop or Tablet PC.