

CODE

PRODUCT DESCRIPTION

Hospital bed Electric - Universal



DA-8-5

**Hospital bed Electric - Universal
with folding sides**


5 FUNCTIONS

High quality cold-rolled steel.
Head and legs power adjustment:
Head :0-75 °, legs :0-40 °
adjustable bed height:
435 mm -850 mm.
Bed head and bed foot tilted :0-20 °
Folding safety rails.
Central Locking swiveling Castor
System, Casters w
Bed head and bed foot removable .
Emergency CPR function.
Bed: 1980 * 900mm,
Bed body: 2090 * 910 mm.
Weight Capacity: 400 lbs(182 kg).
Net Weight: 128KG.
Battery backup



DA-8-1 ICU

Five-function Electric ICU Bed SIZE: L2140×W990×H440-710mm

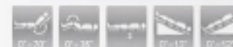



DA-3-5

Hospital bed Electric - Universal with folding sides

3 functions

High quality cold-rolled steel.

Head and legs power adjustment:

Head :0-75 °, legs :0-40 °

Adjustable bed height:

335 mm -750 mm.

Folding safety rails.

Locking swiveling Castor

Bed head and bed foot removable .

Bed: 1980 * 900mm,

Bed body: 2090 * 910 mm.

Weight Capacity: 400 lbs(182 kg).

Net Weight: 108KG.

Battery backup

Box size 2080 x 1050 x 380 mm



DA-3-5

Three-function Electric Bed

SIZE: L2180×W970×H370-680mm



Hospital bed Hydraulic - Universal



S-4

4 - Function hydraulic bed with sides

AIR ASSISTED BACKREST
HYDRAULIC HIGH LOW
TRENDELLENBURG TILTING
REMOVABLE HEAD AND FOOT ENDS
FOLDING BED SIDES

Estimated dimensions 2200 x 950 x 720 mm
WEIGHT ABOUT 120 KG
Bed height with no mattress 500 mm


Four-function Hydraulic Bed



Mesh-wire bed platform is with humanized design, good for ventilation and prevent the mattress from mould.



Trendelenburg and reverse trendelenburg position is controlled by gear mechanic system. Trendelenburg: & Rev. Trend: ±11°

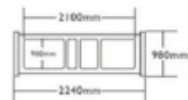

SIZE: L2240 x W980 x H520-810mm

Hospital bed Winding - Universal



A-6-1

3 - Function manual bed
with folding sides

ABS head/foot board

WINDING BACKREST

WINDING KNEE SECTION

FOLDING BED SIDES


REMOVABLE HEAD AND FOOT ENDS

BED HEIGHT WITH NO MATTRESS


400 - 750 MM

Estimated dimensions

2200 x 1100 x 510 mm

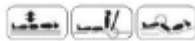

SIZE: L2140 x W940 x H440-700mm

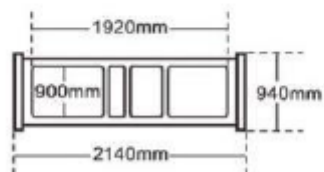



**2-Function full-fowler bed
with folding sides**

ABS head/foot board

WINDING BACKREST

WINDING KNEE SECTION

FOLDING BED SIDES

REMOVABLE HEAD AND FOOT ENDS

BED HEIGHT WITH NO MATTRESS 500 MM

B-16

Estimated dimensions 2100 x 1000 x 250 mm X
2 BEDS

PACKED FLAT

WEIGHT 90 KG ESTIMATED

HEAD AND FOOTENDS X 2 sets (BOX 600 X 250
X 990 mm)

WEIGHT 14.6 KG





B31-1

**1-Function Pneumonia bed with folding sides
With 125 mm castors MA-3D**

Winding Backrest

Folding bed sides

Stainless steel Head and foot-ends

Bed Platform height WITH NO MATTRESS 520
MM

Bed mattress length 1950 mm

Bed mattress width 900 mm

Packing dimensions

2100 mmx 950 mm x 300 mm

2 BEDS per box

WEIGHT 100 KG ESTIMATED





B29

**2-Function Fowler bed with folding sides
With 125 mm castors MA3-D**

Winding operated backrest

Folding bed sides

Stainless steel Head and foot-ends

Bed Platform height WITH NO MATTRESS 520
MM

Bed mattress length 1950 mm

Bed mattress width 900 mm

Packing dimensions

2100 mm x 950 mm x 300 mm

2 BEDS per box

WEIGHT 100 KG ESTIMATED



Standard Hospital beds



BD120

Mobile Hospital bed.

Main frame:

Constructed from 50 x 25 x 1.6mm profile tube.

Mattress seat:

Wire mesh. The mattress seat is reinforced with two lengths of 25 X 5mm flat bar.

Head and foot ends:

Head and foot ends are of steel tube of nominal outside diameter of 31.75mm and wall thickness of 1.6mm. Each end is reinforced with steel tube cross members of nominal outside diameter of 21.34mm and wall thickness of 1.6mm.

Colour and finish:

Epoxy powder coating finish complying with the relevant requirements of SABS1274

Dimensions:

Length 1 980mm, Width: 910mm, Height 150mm PACKED FLAT



BD175

Mobile homecare bed.

Main frame:

Constructed from 50 x 25 x 1.6mm profile tube.

Mattress seat:

Wire mesh. The mattress seat is reinforced with two lengths of 25 X 25mm square tube. 1.6mm wall thickness.

Head and foot ends:

Head and foot ends are of steel tube of nominal outside diameter of 31.75mm and wall thickness of 1.6mm. Each end is reinforced with steel tube cross members of nominal outside diameter of 21.34mm and wall thickness of 1.6mm.

Dimensions:

Length 1 980mm, Width: 910mm, Height 550mm.



BD175

Mobile homecare bed.

Main frame:

Constructed from 50 x 25 x 1.6mm profile tube.

Mattress seat:

Wire mesh. The mattress seat is reinforced with two lengths of 25 X 25mm square tube. 1.6mm wall thickness.

Head and foot ends:

Head and foot ends are of steel tube of nominal outside diameter of 31.75mm and wall thickness of 1.6mm. Each end is reinforced with steel tube cross members of nominal outside diameter of 21.34mm and wall thickness of 1.6mm.

Dimensions:

Length 1 980mm, Width: 910mm, Height 550mm.

BD 175B

SAME BED AS ABOVE BUT WITH A BACKREST

Wire mesh. The mattress seat is reinforced with two lengths of 25 X 5mm flat bar.

Head and foot ends: WITH WOOD INSERTS

Head and foot ends are of steel tube of nominal outside diameter of 31.75mm and wall thickness of 1.6mm. Each end is reinforced with steel tube cross members of nominal outside diameter of 21.34mm and wall thickness of 1.6mm.

Dimensions:

Length 1 980mm, Width: 910mm, Height 150mm PACKED FLAT



BD172

Mobile Hospital bed.

Main frame:

Constructed from 50 x 25 x 1.6mm profile tube.

Mattress seat:

Wire mesh. The mattress seat is reinforced with two lengths of 25 X 25mm square tube. 1.6mm wall thickness.

Head and foot ends:

Head and foot ends are of steel tube of nominal outside diameter of 31.75mm and wall thickness of 1.6mm. Each end is reinforced with steel tube cross members of nominal outside diameter of 21.34mm and wall thickness of 1.6mm.

Colour and finish:

Epoxy powder coating finish complying with the relevant requirements of SABS1274

Dimensions:

Length 1 980mm, Width: 910mm, Height 550mm.



BD 011

Mining Compound and Institutional.

Main frame: Constructed from 45 x 45 x 3mm profile angle.

Mattress seat: Wire mesh or flat rolled expanded metal. The mattress seat is reinforced with two lengths of 25 X 5mm flat bar.

Head and foot ends: Head and foot ends are of steel tube of nominal outside diameter of 31.75mm and wall thickness of 1.6mm. Each end is reinforced with steel tube cross members of nominal outside diameter of 21.34mm and wall thickness of 1.6mm.

Colour and finish: Epoxy powder coating finish

Dimensions:

Length 1 930mm, Width: 910mm, Height 550mm.